


Product
Number **H84102421**
Fits up to 9" neck

Dale® PediPrints™

Dale® 242 PediDucks™ Tracheostomy TubeHolder

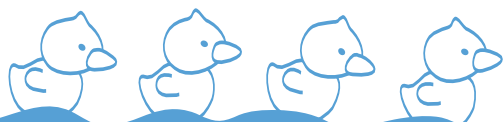

With Moisture Repellent Neckband

LATEX-FREE

DALE®

≡ Patient Care IncorporatedSM



+H84102421.


Product
Number **H84102421**
Fits up to 9" neck

Dale® 242 PediDucks™ Tracheostomy Tube Holder

With Moisture Repellent Neckband





Product Description

Provides secure positioning of the tracheostomy tube. Superior fastener tabs are quickly and easily threaded through the flange of the tracheostomy tube. Two-piece neckband includes moisture repellent lining.

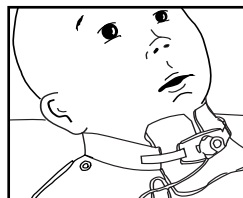

Application Instructions

1. Thread long, narrow fastener tabs through flanges on trach plate.
2. Bring tabs back over flange and adhere to the neckband.
3. Secure fastener tab of short section (A) to the long neckband section (B). Adjust for proper length or compression.
4. Trim excess material (B only) for proper length.

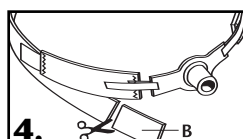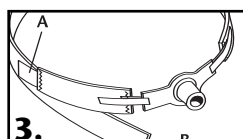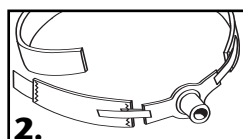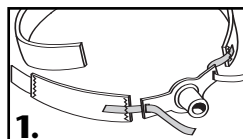

Caution

Use in supervised care settings only. Monitor to prevent unintentional disengagement. Do not trim fastener hook tabs.

For single patient use.

Change holder daily or more frequently as required. Change after patient bathing. Disposable. **Do not wash.**

Ordering Information

- #H8410**2421** Dale® 242 PediDucks™ Tracheostomy Tube Holder
HCPCS Code: A7526 Fits up to 9" neck. Individually packaged. Sold in boxes of 10.
- #H8410**2411** Dale® 241 PediStars™ Tracheostomy Tube Holder
HCPCS Code: A7526 Fits up to 18" neck. Individually packaged. Sold in boxes of 10.
- #H8410**2401** Dale® 240 Blue® Tracheostomy Tube Holder
HCPCS Code: A7526 Fits most. Individually packaged. Sold in boxes of 10.

Dale medical products are available from your preferred hospital/medical supply dealer. **Latex-Free.**

DALE®

DALE MEDICAL PRODUCTS, INC.
7 Cross Street, P.O. Box 1556, Plainville, MA 02762-0556 USA
TEL: 800-343-3980 INT'L: 508-695-9316 FAX: 508-695-6587
www.dalemed.com

CE EU Representative:
MDSS GmbH, Schiffgraben 41
D-30175 Hannover, Germany